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AML1-ETO RQ Kit Manual

**For Real-Time PCR Quantitative Detection of
AML1-ETO Transcripts**

For use with Rotor-Gene or StepOne
Research use only

NG-WI-ASL-31-100
Version 1.1
Winter 2020

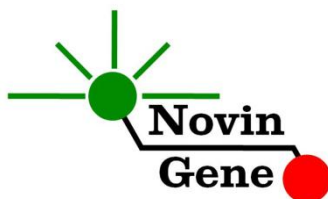
Table of Contents:

1. Introduction	2
2. Kit Contents	3
3. Storage and Stability	3
4. General Precautions	3
5. Additionally Required Materials.....	4
6. Specimen, storage and transport	5
7. RNA isolation	5
8. cDNA synthesis	5
9. PCR Protocol	6
10. Programming of RotorGene	7
11. Programming of StepOne	7
12. Programming of other machines	8
13. Data Analysis: RotorGene.....	8
14. Data Analysis: StepOne.....	11
15. AML-ETO% calculation.....	13
16. Analytical Sensitivity	13

AML1-ETO RQ kit is intended for the detection of AML1-ETO transcripts in peripheral blood and AML1-ETO percentage calculation in patients undergoing therapy. This kit is designed for use with RotorGene or StepOne machines. This kit is for research use only!

Important Note: *This kit doesn't provide reagents for RNA extraction or cDNA synthesis!*

1- Introduction

AML1-ETO is the most common chromosomal aberration in de novo acute myeloid leukemia (AML) patients. This abnormality is resulted from t(8;21)(q22;q22) translocation. The translocation, produces a fusion protein which inhibits myeloid transcriptional factors and as a result, cellular differentiation is blocked. This alteration occurs in approximately 7% of adult and 12% of pediatric with AML.

AML1-ETO RQ provides a ready-to-use system for detection and quantitation of AML1-ETO transcripts as well as calculation of AML1-ETO percentage.

This kit is based on Real-Time PCR technology. In this method application of fluorescent probes allows detection of amplified product. Analysis of fluorescent kinetics also leads to detection of the target sequence in the reaction without requiring post-amplification analysis, reducing the possibility of contamination with the PCR product. This kit is designed to be used with RotorGene or StepOne machines.

2- Kit Contents

The kit contains a manual or a quick guide and following reagents:

Label	Content	Quantity
AML1-ETO Mix*	Master mix for AML1-ETO	480 µl
ABL Mix*	Master mix for ABL	480 µl
AE1	AML1-ETO Standard 1: 100,000 copy/µl	150 µl
AE2	AML1-ETO Standard 2: 10,000 copy/µl	150 µl
AE3	AML1-ETO Standard 3: 1,000 copy/µl	150 µl
AE4	AML1-ETO Standard 4: 100 copy/µl	150 µl
AE5	AML1-ETO Standard 5: 10 copy/µl	150 µl
A1	ABL Standard 1: 100,000 copy/µl	150 µl
A2	ABL Standard 2: 10,000 copy/µl	150 µl
A3	ABL Standard 3: 1,000 copy/µl	150 µl
A4	ABL Standard 4: 100 copy/µl	150 µl
Water	PCR Grade Water	200 µl

* 1, 2 and 4 tubes for 24, 48 and 96 reaction kits.

3- Storage and Stability

The kit components should be shipped and stored at -20°C and are stable until the expiry date mentioned. Avoid repeated freeze-thaws especially for PCR Mix more than few times to prevent reduced sensitivity.

4- General Precautions

In order to prevent false results, always pay attention to the following points:

- Treat all samples as potentially infectious.

- Take utmost care to avoid RNase contamination during RNA extraction and cDNA synthesis.
- Within the pre-PCR work area, assign three separate spaces for: a) Sample storage and extraction; b) Reagent preparation where the Master Mix is aliquoted into tubes; and c) Reaction preparation area for addition of templates to the tubes.
- Always wipe the working surfaces with 70% Ethanol before and after work.
- Thaw kit components on ice completely, mix by flickering followed by a quick spin and store on crushed ice after.
- Do not place 0.2 ml PCR tubes on crushed ice. Use cold blocks instead.

5- Additionally Required Materials

To use this kit, you need the following items:

- Real-Time PCR machine and accessory computer
- Refrigerated microtube centrifuge
- Vortex Mixer
- Dry Block Heater
- Adjustable pipettors and nuclease free filtered tips
- RNA extraction kit
- cDNA synthesis kit
- Nuclease free 1.7ml microtubes and PCR microtubes
- Disposable powder-free gloves
- Cold block

6- Specimen, Storage and Transport

Peripheral blood should be collected in sterile condition in proper and sterile tubes. We recommend EDTA or Citrate as anticoagulant. Whole blood should be shipped and stored at +4°C (stable for 72 hrs).

RNA can directly be extracted from whole blood. Alternately to increase sensitivity, buffy coat can be used. For optimum results, a sample should include 5 to 10 million WBC per 150µl.

Whole blood or buffy coat can be stored at +4°C for three days. Otherwise should be aliquoted and stored at -70°C which is stable for few months.

7- RNA Isolation

RNA isolation can be performed using different kits from various manufacturers. We recommend using:

- TriPure isolation Reagent (Cat. no. 1667157, Roche Applied Science, Germany)
- TRIzol isolation reagent (Cat. no. 15596026, Thermoscientific, USA)
- Isol-RNA isolation reagent (Cat. no. 2302700, 5 prime, Germany)
- Accuzol isolation reagent (Cat. no. K-3090, Bioneer, Korea)

8- cDNA Synthesis

1ug of total RNA is required and should be reverse transcribed to cDNA using random hexamers or specific primers. Different kits are available in the market for this purpose.

These specific primers are available upon request. If the CT observed for ABL using Random Hexamers is above 27, it is recommended to use specific primers for decreasing CT for ABL and AML1-ETO gene and increasing sensitivity of the test.

Dissolve RNA pellet in 5ul of the Primers solution plus 15ul of water, and then add other reagents required for cDNA synthesis.

Dilute prepared cDNA 2.5x with nuclease free water. For example to 20ul of cDNA add 30ul of nuclease free water.

9- PCR Protocol

Thaw the reagents on ice completely, followed by brief mixing and a quick spin.

Each sample should be examined for both of AML1-ETO fusion gene and for ABL control gene. So two set of reactions are required. In AML1-ETO set, consider 1 tube for each sample as well as 6 tubes for the 5 standards (AE1 to AE5) and Negative control or NTC. In ABL set, consider 1 tube for each sample and 5 tubes for the 4 standards (A1 to A4) and Negative control or NTC. Place required number of tubes on cold block.

Pipette 20µl of AML1-ETO Mix to the first series of tubes and 20µl of ABL Mix to each tube of the second group. Continue by adding 5µl of cDNA or standard or control to each tube.

Cap the tubes and visually inspect to make sure all are capped securely. Place tubes in the machine.

Note: Working with StepOne instrument, spin tubes briefly before loading on the block.

Note: If using RotorGene attach the locking ring.

10-Programming RotorGene

Before you start the machine, make sure you have attached the locking ring on the rotor!

Open the CD provided in the kit and double click on “AML1-ETO 0.2” or on “AML1-ETO Strip” depending to tubes used. Click on Start button (Green button on the top menu). On the pop up window click start again and save the run file.

Edit sample names on both AML1-ETO and ABL pages. Remember that tubes containing AML1-ETO Mix should only be named in AML1-ETO page and tubes containing ABL Mix should only be named in ABL page.

Make sure in the “Type” column, all the standards have been defined as “standard” and the relative concentrations have been entered. Patient samples should be defined as “unknown” and no template control as “NTC”, respectively.

11-Programming of StepOne

Open the StepOne software (V 2.*). On the Set Up menu click on Template and select the file on CD provided with the kit. Click on Plate Setup. One negative control, 5 standards for AML1-ETO, 4 standards for ABL and few samples are defined. You may change plate set up using right click options (copy, past, clear). You may also add or remove samples on “Define Targets and Samples” menu. When finished click on Start Run and save the experiment. Instrument will start shortly.

12-Programming Other Machines

If you apply this kit to other Real-Time PCR machines, program it according to the following table:

Step	Temperature and time	Cycles
1	95°C x 10 min	1
2	95°C x 15 sec	45
	60°C x 60 sec	

Fluorescence should be collected at 60°C for FAM and VIC dyes. Both of AML1-ETO Mix and ABL Mix contain ROX with final concentration of 300nM in reaction.

13-Data Analysis: RotorGene

Before analyzing results, make sure in the sample menu all the standards have been defined as “standard” and the relative concentrations have been entered. Patient samples should be defined as “unknown” and Negative control or no template control as “Negative Control” or “NTC”, respectively.

Analyze the data according to manufacturer recommendations. Perform quantitative analysis for both **AML1-ETO (Green channel)** and **ABL (Yellow channel)**. Briefly, click on analysis menu and then, under Quantitation tab double click on “Cycling A. Green”. Close the pop up window and manually set threshold at 0.05. Repeat above for Yellow channel and set the threshold at 0.1. Figures 1 and 2 represent typical graphs for RotorGene machine.

Consider following points when analyzing:

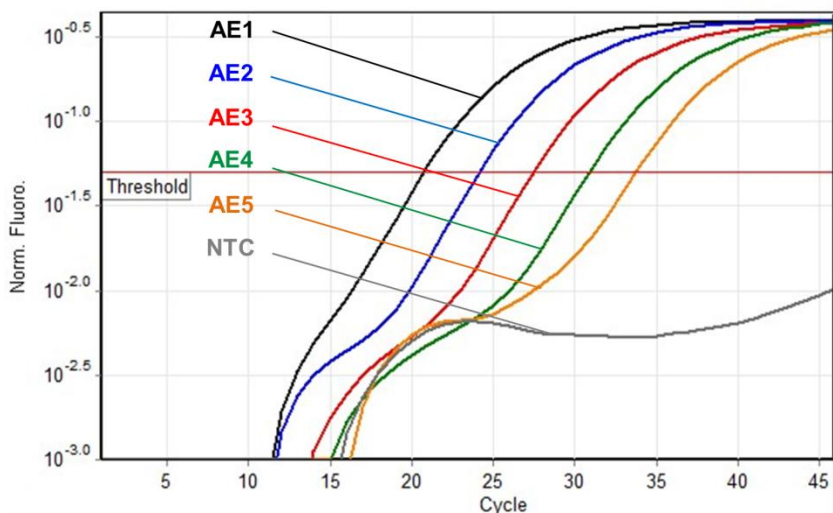
Fig 1. Typical AML1-ETO graph in Green channel for Rotor-Gene

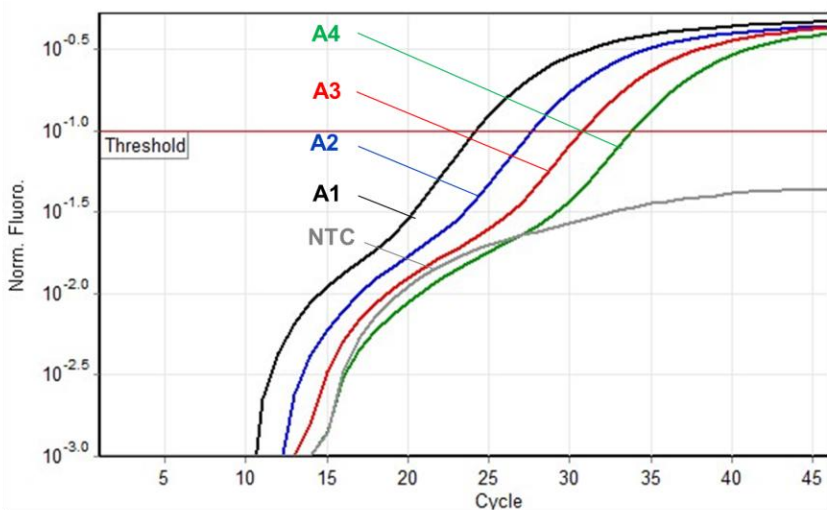
Fig 2. Typical ABL graph in Yellow channel for Rotor-Gene

Note that a sample is considered Positive only if it has a sigmoid graph and log phase, and only then CT is reliable and can be used.

In the absence of sigmoid graph and log phase, sample is considered Negative and CT if present is not reliable.

- A sample is **Positive** if it is positive in both Green/AML1-ETO and Yellow/ABL channels with sigmoid graphs and CT of 20-40 for Green and 20-30 for Yellow.
- A sample is **Negative** if is negative in Green/AML1-ETO channel while it is positive in Yellow/ABL channel with a sigmoid graph and CT of 20-27.
- Results are **inconclusive** and the test should be repeated if a sample is negative in both of Green/AML1-ETO and Yellow/ABL channels. Improper extraction or test set up could cause that.
- Results are **inconclusive** and the test should be repeated if a sample is negative in Green/AML1-ETO channel while it is positive in Yellow/ABL channel with sigmoid graph and CTs above 27. Improper extraction or low RNA input could cause that.

*Note: All patient samples should be positive in Yellow/ABL channel with a sigmoid graph and CT of 27 or less. ABL CT greater than 27 usually happens if not enough cells have been extracted or less than 100ug RNA has been used. CTs higher than 27 reduce the sensitivity of the test and may result in **false negative** reports.*

14-Data Analysis: StepOne

Analyze the data according to manufacturer recommendations. Briefly, click on “Analyze” and set the threshold for **AML1-ETO/FAM** on 0.05 and for **ABL/VIC** on 0.1.

Figures 3 and 4 represent typical graphs for StepOne machine.

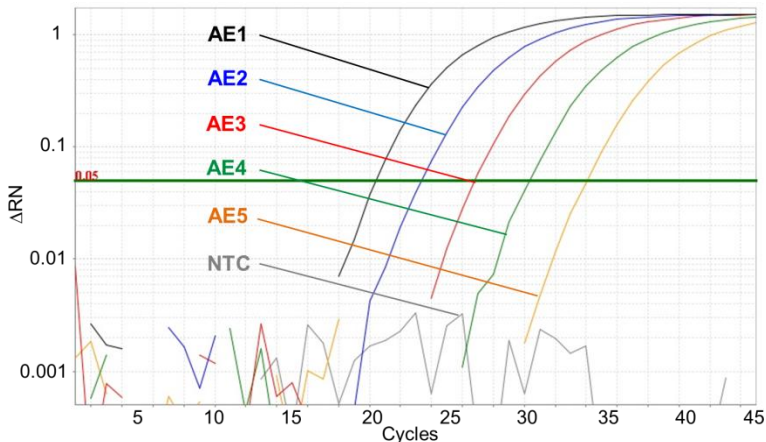
Fig3. Typical AML1-ETO graph in FAM channel for StepOne

Note that a sample is considered Positive only if it has a sigmoid graph and log phase, and only then CT is reliable and can be used.

In the absence of sigmoid graph and log phase, sample is considered Negative and CT if present is not reliable.

Consider following points when analyzing:

- A sample is **Positive** if it is positive in both FAM/AML1-ETO and VIC/ABL channels with sigmoid graphs and CT of 20-40 for Green and 20-30 for Yellow.

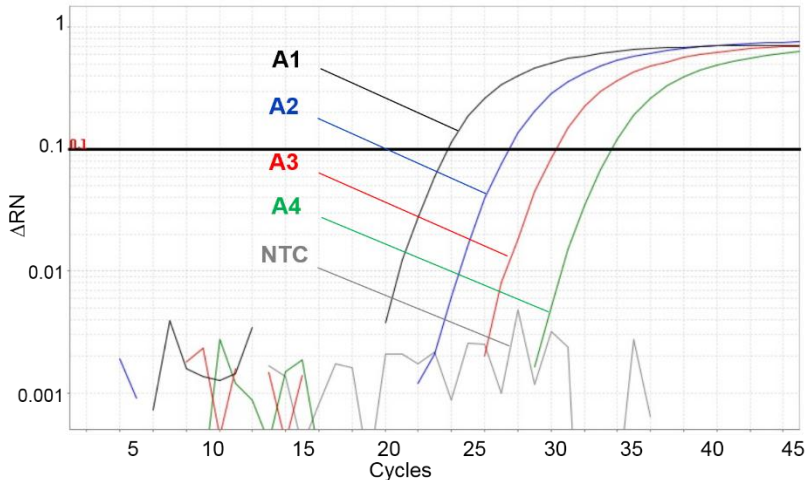
Fig 4. Typical ABL graph in VIC channel for StepOne

- A sample is **Negative** if it is negative in FAM/AML1-ETO channel while it is positive in VIC/ABL channel with a sigmoid graph and CT of 20-27.
- Results are **inconclusive** and the test should be repeated if a sample is negative in both of FAM/AML1-ETO and VIC/ABL channels. Improper extraction or test set up could cause that.
- Results are **inconclusive** and the test should be repeated if a sample is negative in FAM/AML1-ETO channel while it is positive in VIC/ABL channel with sigmoid graph and CT of above 27. Improper extraction or low RNA input could cause that.

*Note: All patient samples should be positive in VIC/ABL channel with a sigmoid graph and CT of 27 or less. ABL CT greater than 27 usually happens if not enough cells have been extracted or less than 100ug RNA has been used. CTs higher than 27 reduce the sensitivity of the test and may result in **false negative** reports.*

15-AML1-ETO% calculation

To assess the response to therapy, AML1-ETO% value for each patient can be calculated. This kit uses NCN method for this purpose (*Beillard E. 2003, Leukemia 17:2474*). In this method AML1-ETO% value is the AML1-ETO transcripts (titer) normalized by the ABL transcripts (titer) and then multiplied by 100%.

Usually titer of ABL gene transcripts is higher than for AML1-ETO fusion gene. Therefore results of above calculation would be a number less than 100%. This value may fall even below 0.001% in optimum cases. However at diagnosis and in case of relapse, AML1-ETO transcripts may surpass ABL gene transcripts. As a result, ratio of AML1-ETO/ABL would be above 100%. Such results have been reported previously too. For example please note table 12 of the paper published by J. Gabert *et al* (2003, *Leukemia 17:2318*) for patients. He reported a range of 110% to 2000% with average of 420% for peripheral blood and range of 29% to 530% with average of 210% for bone marrow samples. Therefore a ratio above 100% is not unexpected and simply means higher transcripts of the fusion gene compared to the control gene.

16-Analytical Sensitivity

Analytical sensitivity of the kit has been determined by examining dilution series of cloned cDNA and was estimated about 2 copy/ul.